



MIYAKO BI PERFECTION

都 美 パーフェクション

The evolution of CaviMAX towards an even flawless beauty

One machine for the perfect

*Body &
Face*



Dual Cavitation &



Moving EMS

Composite beauty instrument for amazing body sensation

LED control panel Hybrid beauty tool !

Multi-use for different needs

Facial
lightening

Body
Shaping

Muscle
care

Facial
care

Lifting
Tightening

Cellulite
care



Ultrasound + Cavitation (Dual Cavitation) + Moving EMS
Triple Function Body Machine


Amazing body sensation beauty instrument equipped with Ultrasound, EMS, and Electric Pulse Method has arrived !!
The oscillation probe, LCD touch screen, High spec total treatment machine

Tones Muscle & Burns Fat Based on science!

Beautiful Body Line from synergistic technology!

Feel the rejuvenating sensation in your body!

Easy to use LCD touch screen!

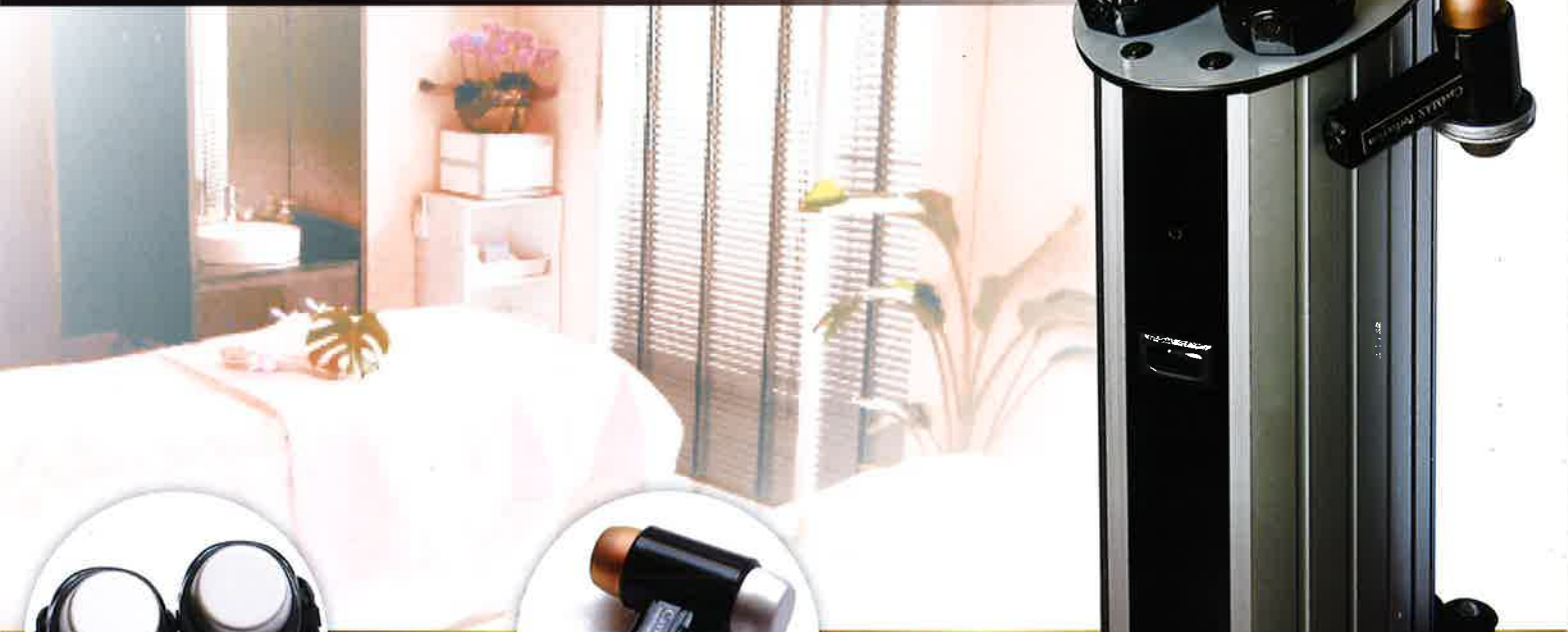


BODY Performance (Slimming Tightening)

- ◆ W Dual cavitation
38 KHz + 1.5 MHz Frequencies
Oscillates simultaneously to stimulate shallow and deep fat deposits
- ◆ W Moving EMS
Complex EMS wave reaches and stimulates deep muscle tissue
- Our O type Probe oscillates with ultrasound technology to create a varied and stable body treatment

Facial performance

- (Facial slimming beautification, Anti-Aging Introduction function)
- ◆ Electro Pulse Method (Penetration)
Injects cosmetic liquid into the deep layers of the skin
 - ◆ 2 MHz Ultra sound
(Facial slimming beautification, Anti-Aging)
2 MHz Ultrasound vibration frequencies are a great solution for slimming the face and reducing the effects of aging
 - ◆ Moving EMS (Tightening and Lift)
EMS will target micro facial muscles, a great solution for tightening and lifting sagging skin due to aging

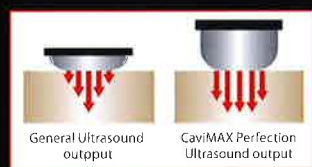


MIYAKO BI PERFECTION
部 大 パーフェクション

All Japanese product

Design Registrations Approved

W Dual cavitation (Ultrasound+cavitation)


W Dual Cavitation is a technology that oscillates 2 separate frequencies. For the shallow fat deposits the beneficial 1.5 MHz frequency is used. For deeper fat deposits, the 38 KHz frequency is used. Both frequencies are combined and generated simultaneously. This combination affects both the shallow and deeper parts of your body adding a warmth that makes the treatment experience comfortable.

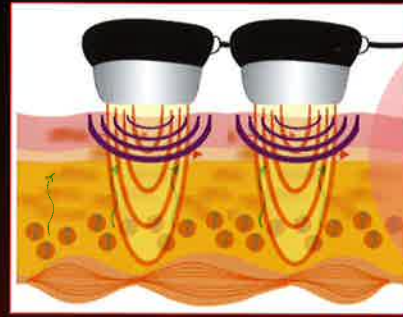

38 KHz cavitation & 1.5 MHz+ Ultrasound frequencies Oscillating simultaneously

W Moving EMS (Low frequency)

By not using a stationary type EMS pad muscular exercises can be achieved from multiple angles. Moving EMS twin head technology will move deep rooted muscle and fascia providing a unique experience in skin care. This simultaneous muscle stimulation is a great treatment for tightening the skin around your hip, stomach, upper arm and back. Our original technology is able to incorporate dual cavitation and EMS oscillation simultaneously using two adjacent probes.

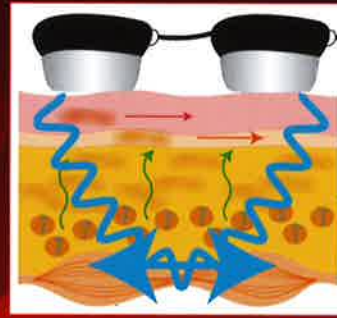

Output from two probes allows for 7 patterns of EMS to move the muscles directly.

CaviMAX Perfection Features

Conventional cavitation can reduce non esterified fatty acids but cannot burn fat. Cavitation alone has always had the fear of fat returning. CaviMAXPerfection combines W cavitation with W moving EMS to further reduce fat through aerobic exercise of the surround muscle. This burns non esterified fat without the fear of fat stubbornly returning.

The ergonomic shape ensures the beauty technicians hands will not get tired during the course of treatment



CaviMAX Perfection Example of use



Waist + Stomach tightening



Hip up +Thigh tightening



Calf tightennng



Back tightening



Upper arm tightening



Bustup

Popular face slimming and lifting effect in tandem with original technology

2 MHz Ultrasound • Moving EMS 【Facial slimming】 (Ultra sound+EMS)

We chose 2 MHz ultrasound frequency as it is known to achieve great results for facial slimming. This frequency also provides comfortable vibrations to the skin.

Moving EMS combined with oscillation simultaneously will cause muscle training and synergize with the effect of the ultrasound treatment. You can use the facial probe for your shoulders and upper arms as well!

Our product will allow the user to choose EMS only, Ultrasound only or both simultaneously.



Facial slimming +Massage



Shoulder + Decollete massage

Electro pulse method (Cosmetic liquid introduction method)

This is the electroporation method also known as mesotherapy. A safe alternative to treatments that require needles. Using our brand new cosmetic liquid application method, we can inject cosmetic gels directly to cells without the need to ionize or polymerize component collagen.



Cosmetic liquid introduction care

※Please clean your skin thoroughly before treatment.
※Safetier VC12、Safetier CLG are sold separately.

Resistant Titanium Head Lasts Over 100 Years!

A genuine titanium head that is safe for your skin!

- Titanium is very gentle and easy on the skin
Like 24k gold, it is a very expensive metal that is not known to cause an allergic reaction. In the medical industry it is used to create prosthetic bones. This genuine titanium probe head will give you better results than products that are simply titanium plated.
- Always hygienic
Titanium is naturally sterilized by light making it always hygienic.
- High safety factor
The structural resistance factor of titanium is very high meaning that the titanium probe is resistance to damage and distortion. This greatly reduces the risk of distortion and buckling of the probes surface, resulting in a product that won't damage important skin.
- Titanium Seamless processing
The T probe and the O probe are seamless and smooth. The probe is built based on human body engineering research and has been designed to fit to all contours of the body. Combined with complete oscillation technology, this probe is able to provide extremely effective skin treatment.



99.6% of titanium is classified as "genuine titanium" by Japanese Industrial Standards: JIS

Exclusive Body Gels' Synergistic Effect

This body gel was created to maximize the cavitation effect. It features 6 important size reducing effects. It also contains many wonderful things for fat control such as: mate leaf extract, colesu root extract, 4 types of caffeine, conjugated linoleic acid, chlorella, fermented white banal pin protein, artichoke extract for blood vessel VP energy, as well as 17 additional slimming components.



※The liquid is pink in color because of the Cyanocobalamin (Vitamin B12) content



- ① T-type Probe (Face use)
- ② O-type Probe x2 (Body use)
- ③ Earth cord for face+Earth titanium plate
- ④ 8 pin cable (Face • Body use)

O-type probe
T-type probe

Ultrasound 1.5 MHz
Ultrasound 2 MHz

Cavitation : 38 KHz
EMS : 3 Hz~280 Hz

EMS : 3 Hz~280 Hz
EP : 1 KHz~25 KHz

Maintenance Guarantee: one year warranty

MIYAKO BI PERFECTION



- Dual Use Face and Body
- 99.6% Titanium head
- O-type probe T-type probe Design
- LCD control display

MIYAKO BI PERFECTION

- Rated voltage AC100V~240V
- Rated frequency 50/60Hz
- Size W280×D360×H780mm
- Electricity consumption 36W
- Weight 13.0kg

Exclusive cosmetics



- ① Exclusive face gel
① Safetier gel cream plus (200g)
- ② Exclusive face gel
② Safetier gel cream plus (980g)
980g Exclusive gel slide holder
- ③ Exclusive body gel
③ Safetier body gel SP (1kg)



- 1, EP Carrier jel (100g)
- 2, EP essence Ar (20ml) (Arbutin)
- 3, EP essence Vc (20ml) (Triple Vitamin)
- 4, EP essence Ac (20ml) (Argireline)
- 5, EP essence Pc (20ml) (Placenta)

Cosmetic liquid (Electro Pulse Method)
④ EP essence series

Safety handling instructions (Please read the operating manual and all other safety instructions before using the equipment)



Warning

You should not use this product if you have following medical conditions

- Heart disease ● Acute disease ● Contagious disease ● Malignant tumor ● Those of anemia ● Febrile disease ● Blood pressure issues ● Epilepsy ● Expectant and nursing mothers ● Skin conditions ● If you feel any abnormality

● If you develop any abnormal symptoms from using this product. Please consult your doctor before using this product

Manufacture Agent

Miyako Chemical Co., Ltd.
19 Hiramachi Saiin Ukyo-ku Kyoto City Japan 615-0022
Ph:+81 75 312 7247 Fax:+81 75 322 7002

Sales Agent

AVANT BEAUTY Co., Ltd.
KOWA Bld.4F 5-1-12 Shiba Minato-ku Tokyo Japan 108-0014
Free Dial:+81 120 125 222 Fax:+81 3 5484 7521